

Product Information - Autoclaveable HA

LIFE SCIENCES

Autoclaveable HA Pressure Transmitter

Introduction

Anderson-Negele's Autoclaveable HA Sanitary Pressure Transmitter is specifically designed for critical measurements in applications where equipment is autoclaved. Developed from the highly stable "Mini" HA Pressure Transmitter, the Autoclavable HA takes the versatility of a small form factor transmitter to the next level by withstanding the stress of the autoclave environment. Its one-piece stainless steel package design incorporates a transducer and electronic circuitry to convert pressure and/or vacuum to a proportional 4-20 mA signal.

The product can be ordered with 1-1/2" Tri-Clamp® connection or our CPM fitting that eliminates hold-up volumes in process lines. All units are supplied factory calibrated to standard or custom ranges. For field maintenance, non-interactive zero, and span adjustments, as well as easily accessible test points that greatly simplify calibration and validation. 316 "L" stainless steel wetted parts are electro-polished to an Ra=8. Finally, since it's an Anderson-Negele "Life Sciences Series" product, every transmitter includes calibration and material certificates and permanent tag at no charge.

Complete specifications and ordering information are available on the reverse. For more information please visit our Web Site at www.anderson-negele.com, or contact your local Authorized Anderson-Negele Distributor.



Features

- · Compact low profile design
- · Output remains stable through repeated autoclave cylces
- · Rated NEMA 4X & 6P, IP67 & 69k

LIFE SCIENCES

Specifications

Excitation: 10-40 VDC (Absolute), 24 VDC

Nominal regulated or unregulated

Output: 4-20 mA DC, 2 wire with non-

interrupting circuit verification test points

Loop Resistance: 0-700 ohms at 24 VDC

Wiring Connection: 5pin M12 Quick Disconnect Receptacle Recommended Cable: 22-24 AWG, foil shielded, and PVC coated.

(3/16 - 1/4 OD insulation)

Accuracy: ± 0.4% of full scale
Repeatability: ± 0.2% of full scale
Hysteresis: ± 0.10% of full scale
Linearity: ±0.10% of full scale

Stability: ±0.15% of calibrated range for 6 months or

30 Autoclave cycles

Effect of Temperature

Change: ± 0.1 psig/5.5°C (10°F) typical

Max. Autoclave

Temperature: 124°C (255°F) for 1 hour

Process Temperature

Limits: -1°C to 149°C (30°F to 300°F)

(Horizontal Mount recommended over

135°C (275°F))

Ambient Temperature

Limits: -18°C to 49°C (0°F to 120°F)
Storage Temperature: -40°C to 65°C (-40°F to 149°F)

Over-Range Rating: 2 times base range

Response Time: 200 uSec

Wetted Parts: 316L stainless steel electropolished

(Ra max. = 8 microinches, .2 microns)

Housing Material: 304 stainless steel

Actuating Fill: 100% mineral oil. Meets FDA requirements

(21 CFR, 172.878 and 178.3620(a))

Zero and Span

Adjustments: ± 10% of range
Mounting: Direct connection

Electrical Ratings: Intrinsically safe for use in Class 1, Div. 1,

Groups A-D; CE

Compliant: (EMC Industrial Locations)


Standards: Designed and manufactured to sound

engineering practices in accordance with

Article 3.3 of the PED 97/23/EC

NEMA 4X & 6P IP67 & 69k CSA B51-03

Order Information

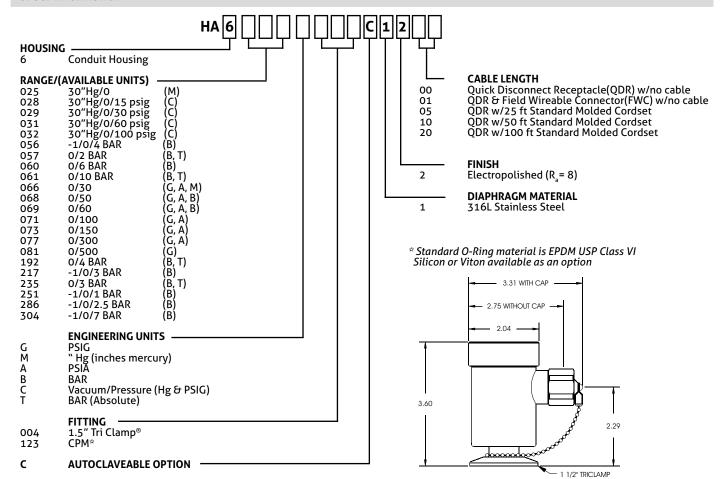

05049 / 2.0 / 2015-03-20 / PW / NA

ANDERSON INSTRUMENT COMPANY 156 Auriesville Road Fultonville, NY 12072 Phone 800-833-0081 Fax 518-922-8997 info@anderson-negele.com Tech. Support: techservice@anderson-negele.com Phone 800-833-0081

SHOWN WITH AUTOCLAVE CAP